


Wound Drainage System (100% Silicon)

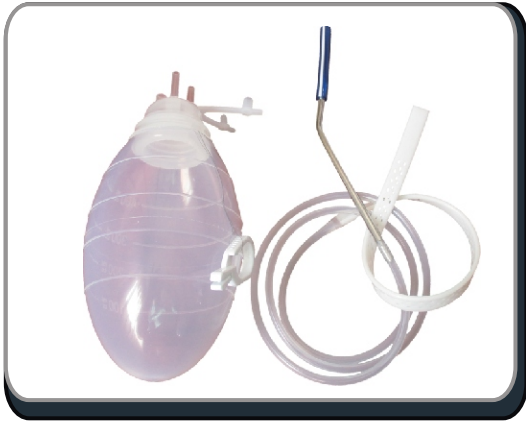

Code : SR-08

- ◆ It consist of silicon make bellow unit. Is manufactured from 100% silicone & is latex free. In-built non-return valve prevents reflex of fluid.
- ◆ Bellow provided with a convenient clip for secure hanging drain.
- ◆ 100% silicone & latex-free , 20cm long, perforated, 7mm x 3mm flat drain, that enables easy placement & prevents accidental sliding out.
- ◆ Radio opaque feature of drain facilitates accurate placement.
- ◆ Flat Drain attached to a 60cm long 100% silicone tubing for attachment to the Bellow unit. Flat drain causes minimum trauma during removal.
- ◆ Unit Consist of the following items :
Silicon Bellow, Silicon perforated catheter, Stainless steel trocar
- ◆ Available Sizes : 100ml , 200ml , 400ml & 800ml.
- ◆ Individually packed in medical grade paper plastic pouch.